

Bench-scale experiments for the development of a unified loop-mediated isothermal amplification (LAMP) assay for the *in vitro* diagnosis of *Leishmania* species' promastigotes

M. KARANI¹, I. SOTIRIADOU^{1,2}, J. PLUTZER³ AND P. KARANIS^{1*}

¹ Medical and Molecular Parasitology Laboratory, University of Cologne, Medical School, Centre of Anatomy, Institute II, Cologne, Germany

² Centre of Dental Medicine, Polyclinic of Operative Dentistry and Periodontology, University of Cologne, Cologne, Germany

³ National Institute of Environmental Health, Department of Water Hygiene, Gyáli ut 2-6, Budapest, Hungary

Received 24 March 2013; Final revision 1 September 2013; Accepted 25 September 2013;
first published online 29 October 2013

SUMMARY

We developed, in bench-scale experiments, a unified loop-mediated isothermal amplification (LAMP) assay for the detection of cutaneous, mucocutaneous and visceral leishmaniasis using DNA of cultivated promastigotes. Two primer sets for the LAMP assay were designed based on the 18S rRNA gene, and their sensitivity and specificity were tested and compared. Both of them were specific for *Leishmania* as the DNA of all ten *Leishmania* species tested was amplified, whereas the DNA of other parasites, including that of *Trypanosoma*, was not. The detection limit for primer set 1 ranged between 30 pg and 3·6 fg, depending on which *Leishmania* species tested. Primer set 2 showed high sensitivity, but was less sensitive than primer set 1. Our findings lead to the conclusion that the LAMP assay with primer set 1 is a promising and effective assay for the successful detection of a wide range of *Leishmania* infections using only a unified multiplex LAMP test.

Key words: LAMP, leishmaniasis, molecular biology, parasites.

INTRODUCTION

Leishmaniasis refers to a diverse group of diseases caused by *Leishmania* spp., protozoan parasites of the family of Trypanosomatidae. Leishmaniasis is geographically widespread because different *Leishmania* spp. are capable of adapting to different environmental circumstances. *Leishmania* spp. have established themselves in urban and rural settings, in deserts and forests and are transmitted by the bite of an infected sand fly. Humans exhibit visceral, cutaneous and mucocutaneous forms of infection [1]. The

variability of the disease depends on the host's immune response and on the particular species of *Leishmania* [2–4].

Leishmaniasis is endemic in 98 countries worldwide, and more than 350 million people are exposed to infection [1]. It has been recently estimated that 0·5 million new cases of visceral leishmaniasis occur each year, with 50000 cases leading to death, and 1·5 million cutaneous leishmaniasis cases occur annually [1].

The diagnosis of leishmaniasis in humans is usually made by serological tests, *in vitro* cultivation and detection of the parasites in tissue aspirates [5, 6]. In the last two decades, Polymerase chain reaction (PCR) assays have been used for more sensitive and accurate diagnosis of these diseases [7, 8]. However, all

* Author for correspondence: Prof. Dr. P. Karanis, University of Cologne, Medical School, Centre of Anatomy, Institute II, Cologne, Germany.
(Email: panagiotis.karanis@uk-koeln.de) (P. Karanis)

of these methodologies are time-consuming and require special, expensive equipment and a well-trained and experienced staff [9]. These circumstances are lacking in developing countries, particularly in the areas where the disease is most widespread. Furthermore, in subclinical cases or when the number of parasites is low, such as in patients with visceral leishmaniasis, microscopy may be unsuccessful in detecting the disease. Taking these facts into consideration, there is a great demand for a technique that will enable a definite, simple, rapid and cost-effective diagnosis of leishmaniasis.

In recent years, a new molecular-based method, i.e. loop-mediated isothermal amplification (LAMP), has been successfully applied for the detection of other pathogens. LAMP promises to be a rapid, specific, sensitive and simple assay for pathogen detection and identification. LAMP can amplify DNA or RNA under isothermal conditions in 1 h [9] and its popularity has rapidly grown for different applications and areas [10]. LAMP has already been used for the diagnosis of several parasitic diseases, including malaria [11], cryptosporidiosis [12, 13], toxoplasmosis [14], trypanosomiasis [15] schistosomiasis [16].

The aim of the present study was to develop and evaluate a unified LAMP protocol in bench-scale experiments *in vitro* and to utilize it for the detection of a wide range of *Leishmania* spp. that cause cutaneous, mucocutaneous and visceral leishmaniasis. Two LAMP primer pairs for the amplification of *Leishmania* DNA were designed to provide improved coverage of globally diverse species together with high specificity as no other protozoan DNA was amplified.

MATERIALS AND METHODS

Origin of *Leishmania* DNA

The development of LAMP assays required DNA from *Leishmania* spp. Seven reference strains were cultivated under sterile conditions and incubated at 26 °C. The *Leishmania* spp. that were used for cultivation were *L. donovani* (zymodeme MON-2, MHOM/IN/1996/THAK35), *L. major* LV39 (zymodeme MON-4, MRHO/SU/59/P), *L. infantum* L4 (zymodeme MON-1) (obtained from the Pasteur Institute, Athens, Greece) and *L. major* 230/119 (zymodeme MON-4), *L. major* 159/06, *L. tropica* 208/48 (zymodeme MON-60), *L. tropica* 101/65 (obtained from the Institute of Medical Parasitology, Koblenz, Germany).

Three types of media were prepared and used for the cultivation of the individual strains. *L. major* LV39 and *L. donovani* MHOM/IN/1996/THAK35 were cultivated in RPMI-1640 medium [17], whereas *L. infantum* L4 was cultivated in DMEM. Leibovitz's medium was used for the cultures of *L. major* 230/119, *L. major* 159/06, *L. tropica* 208/48 and *L. tropica* 101/65 [18]. *Leishmania* subcultures were generated by inoculating tubes containing fresh medium with parasites once a week under sterile conditions.

The DNA of the cultivated parasites was extracted by the phenol-chloroform extraction method [18] and frozen at –20 °C until use.

DNA from *L. guyanensis* (referred to by the WHO as MHOM/GF/79/LEM85), *L. braziliensis* (referred to by the WHO as MHOM/BR/00/LTB300), *L. donovani* DON-43, *Trypanosoma evansi*, *T. brucei*, *Cryptosporidium parvum*, *Toxoplasma gondii*, *Giardia lamblia* and *Entamoeba histolytica* were also used for conducting specificity tests.

Primer design

Two sets of four primers were designed to be specific for *Leishmania* spp. based on the 18S rRNA gene (Table 1, Figs 1 and 2). The sequence of the *L. major* rRNA-encoding gene on chromosome 27, from nt 989660 to 991863 was used. The primers were designed to detect *L. donovani*, *L. major*, *L. infantum*, *L. braziliensis*, *L. guyanensis*, *L. tropica*, *L. mexicana*, *L. amazonensis*, *L. chagasi*, *L. sp. siamensis*, *L. naiffi*, *L. lainsoni* and *L. panamensis*.

LAMP

The components of the reaction mixture were 1 µl primer mix (5 pmol F3, 5 pmol B3, 40 pmol FIP, 40 pmol BIP), 1 µl (8 U) *Bst* DNA polymerase, 2 µl DNA sample, 12.5 µl of 2x LAMP buffer containing 1.4 mM each of dNTP, 0.8 M betaine, 20 mM Tris-HCl (pH 8.8), 10 mM KCl, 10 mM (NH₄)₂SO₄, 8 mM MgSO₄ and 0.1% Tween 20. The mixture was tested with both primer sets, and the reactions were performed at temperatures ranging from 58 °C to 70 °C and for incubation times ranging from 50 to 120 min to determine the optimal amplification conditions. Visualization of LAMP products was performed by gel electrophoresis using 1.6% agarose gels stained with ethidium bromide (1 µg/µl) and run for 30 min at 100 V in an electrophoresis chamber. When the electrophoresis was complete, the gel was carefully

Table 1. Nucleotide sequences of LAMP primers (sets 1 and 2) designed for detection of *Leishmania* promastigotes

Assay	Primer type	Primer length	Primer Sequence	Amplicon size	Target gene
LAMP set 1	FIP (F1c-F2)	39	5'-CACCGGACGCACTAAACCCC-GTGCTTTCCCATCGCAACT-3'	198	FR796423 <i>Leishmania</i> major strain Friedlin complete genome, chromosome 27 from nt 990604–990773
	BIP (B1-B2c)	39	5'-TCCGGCCGTAACGCCTTTTC-CAGTACGTTCTCCCCCGAA-3'		
	F3	20	5'-ATTATGGAGCTGTGCGACAA-3'		
	B3	20	5'-AGAATTTACCTCTGACGCC-3'		
LAMP set 2	FIP (F1c-F2)	41	5'-TGGGTAAAACCCGCCGATGAG-GCTGCCCAGTAGAATTCAGA-3'	190	FR796423 <i>Leishmania</i> major strain Friedlin complete genome, chromosome 27 from nt 990305 to 990492
	BIP (B1-B2c)	42	5'-GGTCGAATTCTTCTCTGCGGGA-CATCACAGACCTGCTGTTGC-3'		
	F3	19	5'-CGTCAACGGACGAGATCCA-3'		
	B3	19	5'-CGCCAGAACATTGAGGAG-3'		



Fig. 1. LAMP primer set 1, targeting the *Leishmania major* rRNA-encoding gene on chromosome 27. The partial sequence of *L. major* and the location of four primers, FIP (F1c-F2), BIP (B1c-B2), F3, and B3 are shown. The arrows indicate the direction of the extension. The numbers on the left indicate the nucleotide position.

transferred to a gel documentation system and visualized by UV transillumination.

Nested PCR amplification

For the PCR detection of *Leishmania*, the gene selected was the internal transcribed spacer-1 gene (ITS-1), which is part of the ribosomal RNA gene of *Leishmania*. The primers used were LITSRn (5'-CTG GAT CAT TTT CCG ATG-3') and L5.8s (5'-TGA TAC CAC TTA TCG CAC TT-3') [18]. The reaction

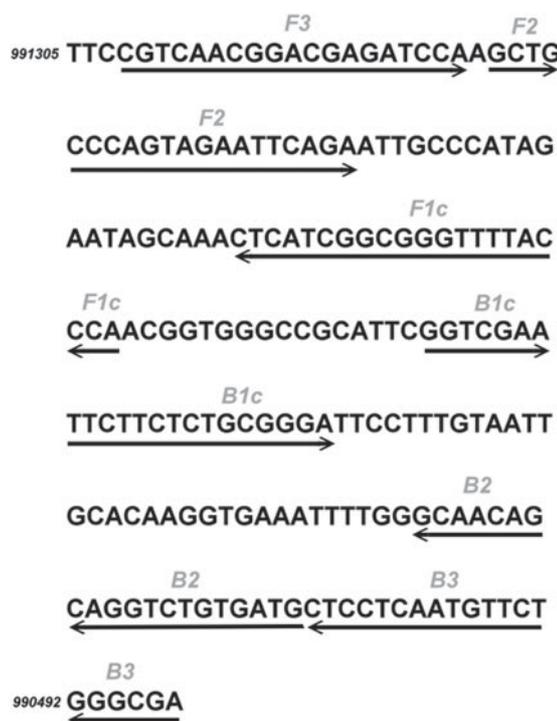


Fig. 2. LAMP primer set 2, targeting the *Leishmania major* rRNA-encoding gene on chromosome 27. The partial sequence of *L. major* and the location of four primers, FIP (F1c-F2), BIP (B1c-B2), F3, and B3 are shown. The arrows indicate the direction of the extension. The numbers on the left indicate the nucleotide position.

mixture contained 5 μ l of 10x PCR buffer (Qiagen, Germany), 2.5 mM each of dNTP (Fermentas, Canada), 10 pmol of each primer, 1 U HotStarTaq polymerase (Qiagen), 2.5% DMSO, distilled water and 5 μ l DNA template. The conditions of the amplification in the thermal cycler were 95 $^{\circ}$ C for 15 min, then 34 cycles of the following steps: denaturation at

95 °C for 20 s, annealing at 53 °C for 30 s and extension at 72 °C for 60 s. After completion of 34 cycles, the reaction was subjected to final extension of 72 °C for 6 min. The second PCR was performed under the same amplification conditions and with the same reaction mix as the first PCR, except that the DNA template was 1 µl PCR product from the first amplification. The PCR products were electrophoresed as described above.

Sensitivity

The sensitivity of the PCR assay was tested using samples containing the target DNA of *L. major* 230/119 with concentrations of DNA ranging from the DNA of 1×10^5 promastigotes/5 µl DNA template to the DNA of 1 promastigote/5 µl. Promastigotes were counted by haemocytometer prior to DNA extraction. A negative control of a sample containing distilled water instead of DNA was always included.

Specificity

To ensure that DNA from several *Leishmania* spp. and strains could be amplified, the specificity of the PCR assay was assessed using cultured reference strains of *L. donovani* MHOM/IN/1996/THAK35, *L. infantum* L4, *L. major* 159/06, and *L. major* LV39. Additionally, the DNA of various other parasites, including *T. evansi*, *T. brucei*, *Toxoplasma gondii*, *C. parvum*, *G. lamblia* and *E. histolytica*, was assayed to test the specificity of the PCR assay.

LAMP primer working conditions

The LAMP product was detected by agarose gel electrophoresis and exhibited a specific pattern of multiple bands. After examining all the cultivated reference strains, the optimal conditions for amplification were found to be 65 °C and 85 min for primer set 1 and 63 °C and 85 min for primer set 2. At the end of the reaction, the temperature was increased to 80 °C for 5 min to terminate the reaction.

RESULTS

Testing LAMP primer set 1

The DNA from ten different *Leishmania* isolates as well as from other parasites was examined to test the specificity of primer set 1. All samples containing DNA from the *Leishmania* spp. (*L. donovani* THAK35, *L. infantum* L4, *L. major* LV39,

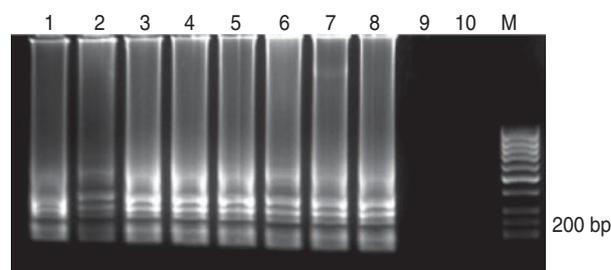


Fig. 3. Specificity of the LAMP assay using primer set 1 on different *Leishmania* isolates at 65 °C for 85 min. Lane 1, *Leishmania major* 159/06; 2, *L. major* 230/119; 3, *L. infantum*; 4, *L. tropica* 208/48; 5, *L. tropica* 101/65; 6, *L. donovani* THAK35; 7, *L. guyanensis*; 8, *L. braziliensis*; 9, *Trypanosoma evansi*; 10, negative control (distilled water); M, 100 bp DNA ladder.

L. donovani DON-43, *L. tropica* 208/48, *L. tropica* 101/65, *L. guyanensis*, *L. braziliensis*, *L. major* 159/06, *L. major* 230/119) were detected by primer set 1 (Fig. 3). The DNA from other protozoa (*T. evansi*, *T. brucei*, *C. parvum*, *G. lamblia*, *E. histolytica*) was not amplified, ensuring that differential diagnosis could be made with these primers.

The sensitivity of the primers was estimated with serially diluted DNA samples from various *Leishmania* spp. Under the conditions of 65 °C for 85 min, the detection limit was variable for each *Leishmania* strain, ranging from 30 pg to 3.6 fg. The detection limits are shown in Table 2. The highest sensitivity was for the sample of *L. tropica* 101/65, for which the detection limit was found to be 3.6 fg, followed by *L. major* 230/119 at 12.6 fg and *L. infantum* L4 at 46 fg. The sensitivity for the sample of *L. donovani* THAK35 was determined to be the lowest at 30 pg, followed by the samples of *L. major* and *L. tropica*.

Testing LAMP primer set 2

The DNA of the same *Leishmania* isolates and heterologous DNA from other protozoan parasites were also tested using primer set 2 to determine the specificity of the primer set. DNA from all *Leishmania* isolates was amplified, whereas the DNA of the other parasites failed to provide a positive result. Therefore, primer set 2 is also specific for the detection of *Leishmania* spp.

Similar to the previous trials with primer set 1, after examining different amplification times and temperatures using all 10 strains, the optimal conditions for amplification were found to be 63 °C and 85 min. To examine the sensitivity of LAMP primer set 2,

Table 2. Sensitivity of LAMP primer set 1 using different *Leishmania* strains under optimum amplification for set 1 (65 °C for 85 min) and for set 2 (63 °C for 85 min)

<i>Leishmania</i> species	<i>L. donovani</i> THAK35	<i>L. infantum</i> L4	<i>L. major</i> LV39	<i>L. tropica</i> 101/65	<i>L. tropica</i> 208/48	<i>L. major</i> 230/119	<i>L. major</i> 159/06
Primer set 1							
Detection limit	30 pg	46 fg	160 fg	3.6 fg	80 fg	12.6 fg	62 fg
Primer set 2							
Detection limit	30 pg	46.2 pg	1.6 fg	3.6 pg	n.d.	n.d.	n.d.

n.d., Not done

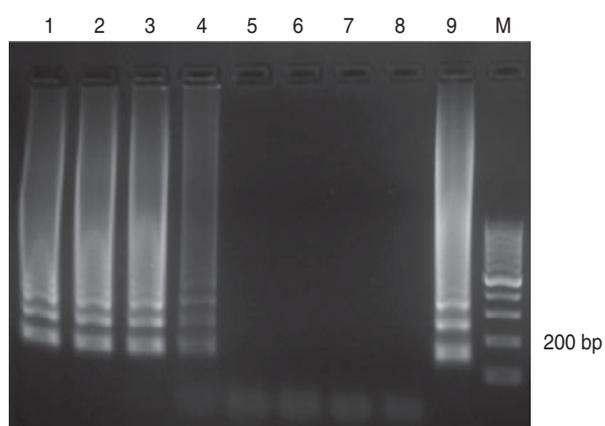


Fig. 4. Sensitivity of the LAMP assay using primer set 2 on serially diluted DNA of *Leishmania tropica* 101/65 at 63 °C for 85 min. Lane 1, 3.6 ng; 2, 360 pg; 3, 36 pg; 4, 3.6 pg; 5, 360 fg; 6, 36 fg; 7, 3.6 fg; 8, negative control; 9, positive control (*L. infantum*); M, 100 bp DNA ladder.

the lowest concentrations of the amplified strain of each subtype of *Leishmania* that grew best during the previous trials were tested, namely *L. donovani* THAK35, *L. infantum* L4, *L. major* LV39 and *L. tropica* 101/65 (Fig. 4) at 63 °C for 85 min. The results, shown in Table 2, indicate that primer set 2 is sensitive for the *Leishmania* spp. tested. The sensitivity for *L. donovani* THAK35 was 30 pg, and the sensitivity of *L. infantum* L4 and *L. tropica* 101/65 was 46.2 pg and to 3.6 pg, respectively. The only sample that was detected with greater sensitivity using this primer set was *L. major* LV39, with a detection limit of 1.6 fg.

Results for PCR

PCR was performed, and the size of the PCR product was found to be 300–350 bp. This product represents the amplified part of the ribosomal ITS-1 DNA fragment. Using PCR, we successfully amplified DNA extracted from $1 \times 10^5/5 \mu\text{l}$ to 1×10^1 parasites/ $5 \mu\text{l}$ of *L. major* DNA template. The results were visualized

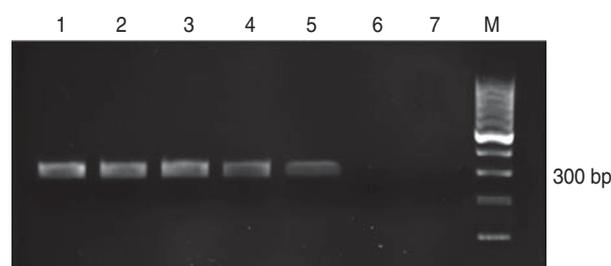


Fig. 5. PCR assay on serially diluted *Leishmania major* promastigotes. Lane 1, 105 promastigotes; 2, 104 promastigotes; 3, 103 promastigotes; 4, 102 promastigotes; 5, 10 promastigotes; 6, 1 promastigote; 7, negative control; M, 100 bp DNA ladder.

by gel electrophoresis and are shown in Figure 5. A DNA band of 300–350 bp can be seen in the first five lanes, whereas there is no amplification seen in lanes 6 and 7, which represent one parasite and the negative control, respectively. Considering that 1 fg is equivalent to ~ 0.1 parasite [19], the detection limit of nested PCR was ~ 100 fg.

All of the examined DNA samples from different *Leishmania* spp. were successfully amplified by PCR assay. The DNA of the other tested organisms was not amplified, indicating that this assay was specific only for species of the *Leishmania* genus.

DISCUSSION

LAMP is a novel molecular method that may be used as an alternative method for the diagnosis of leishmaniasis [20]. We developed and evaluated LAMP assays with two primer sets for the detection of a wide range of *Leishmania* spp. using DNA extracted from cultivated *Leishmania* promastigotes. Both primer sets were specifically designed to detect *Leishmania* spp. using only one test. These results indicate that different samples of *L. major* and *L. tropica* have variable sensitivity. However, it can be

concluded that these primers are specific only for species of the *Leishmania* genus because they produced no false positive results for heterologous species in the LAMP assay. The DNA of 10 *Leishmania* spp. was amplified by both LAMP assays, whereas that of other parasites, including *T. evansi* and *T. brucei*, was not, indicating the specificity of both primer sets for the *Leishmania* genus. The control of the reaction with *Trypanosoma* DNA is of great importance because of the phylogenetic similarity between *Leishmania* spp. and *Trypanosoma* spp.

For primer set 1, the detection limit was different for each species, ranging from 30 pg to 3.6 fg. The *Leishmania* parasitaemia threshold is in the range of 1–8 parasites/ml [19]. The sensitivity was at the fg level, except for *L. donovani* THAK35. A possible explanation could be the presence of genetic heterogeneity in the subtypes. Primer set 2 also exhibited high sensitivity. However, its sensitivity was lower than that of primer set 1, with the exception of *L. major* LV39, for which the detection limit was 1.6 fg (a 100-fold increased sensitivity compared to primer set 1). In DNA samples containing *L. infantum* L4 and *L. tropica*, LAMP primer set 2 exhibited a 100-fold decrease in sensitivity compared to primer set 1. For the isolate *L. donovani* THAK35, the results were identical using both primer sets when an equal amount of DNA was amplified. Therefore, this primer pair could also be suitable for the detection of *Leishmania* DNA. The detection limit of the PCR assay was 1×10^1 parasites of *L. major*. As mentioned above, considering that 1 fg is equivalent to ~ 0.1 parasite [21], the detection limit of PCR was ~ 100 fg. Therefore, LAMP primer set 1 was more efficient than primer set 2.

The high sensitivity and specificity of LAMP primer set 1 indicates that it is an effective method for the successful detection of all pathogenic species of *Leishmania*. However, further trials must be undertaken, as our trials were performed with DNA from cultured promastigotes from the *in vitro* system. An evaluation with tissue samples is necessary so that this technology can be used to diagnose human and animal leishmaniasis, as LAMP-based point-of-care diagnostics should support the diagnosis of clinical leishmaniasis cases.

A LAMP assay for the amplification of *L. donovani* has been reported by Takagi *et al.* [22] with a detection limit of 1 fg. The sensitivity appears to be high, but this assay was tested only with one *Leishmania* spp. In another study, a *Leishmania* assay was

developed for the detection of several species in one test [23]. The sensitivity was between 10 and 100 parasites/ml, using an RT–LAMP assay with loop primers [23]. A disadvantage of the pan-RT–LAMP assay was the observation of false positives from *T. brucei* and *T. cruzi* due to the similarity of the 18S region between members of the Kinetoplastida class. The authors tried to underestimate the degree of cross-reactivity based on the different symptomatology of Chagas disease and leishmaniasis. They recommended the design of species-specific LAMP reactions or degenerate primer sets [23].

In the present paper, we describe a protocol for LAMP that is capable of amplifying a wide range of *Leishmania* spp. in a single test. PCR has been found to be a very good diagnostic method, and many studies of leishmaniasis are based on PCR. However, LAMP only requires basic equipment (no more than a water bath, when first evaluated and standardized) as the temperature is constant; it is therefore easier to perform in countries where leishmaniasis is endemic and costs must be minimized [23]. LAMP is a cheap test as any developed test has its own costs. Once the method is considered as developed and evaluated, the cost for the application of the method is significantly lower. We compared LAMP with a well optimized and sensitive PCR protocol which has been previously published [18]. Large hospitals that possess a real-time turbidimeter would have the benefit not only of a fast and sensitive diagnosis in complicated or subclinical cases but also the opportunity to use LAMP as a quantitative method. In addition, a turbidimeter could be used for the diagnosis of a wide range of diseases. In the future, we expect that LAMP assays will accelerate investigations for the detection of *Leishmania* infections in humans and reservoir animals, including epidemiological/entomological investigations in sand flies, as material from infected sand flies can be used for LAMP assays. Further research is required for a definite assessment of the application of LAMP as a supplement or routine diagnostic method for the detection of *Leishmania* infection in hosts and vectors [24].

In conclusion, a new, rapid, simple and sensitive LAMP method for the diagnosis of *Leishmania* infection has been developed. We anticipate that this unified LAMP assay will contribute significantly to the diagnosis of the disease and can improve the convenient treatment and epidemiological studies of leishmaniasis in any setting worldwide, even if it does not differentiate between species of *Leishmania*. It is

expected that this protocol will find applications in various countries using material from different sources and confirm the current assay as a useful and effective tool in the routine diagnosis of leishmaniasis.

ACKNOWLEDGEMENTS

We thank Dr Eleni Dotsika (Pasteur Institute Athens, Greece), Dr Patrick Scheid (Institute Medical Parasitology, Koblenz, Germany) and Dr Gabriele Schönian (Institute for Mikrobiologie und Hygiene, Charité Berlin, Germany) for supporting this study by supplying us with *Leishmania* parasites and DNA.

DECLARATION OF INTEREST

None.

REFERENCES

1. **World Health Organization.** Control of the leishmaniasis: report of a meeting of the WHO Expert Committee on the control of Leishmaniases. Geneva, 2010.
2. **Anstead G, et al.** Malnutrition alters the innate immune response and increases early visceralization following *Leishmania donovani* infection. *Infection and Immunity* 2001; **69**: 4709–4718.
3. **Olivier M, Gregory DJ, Forget G.** Subversion mechanisms by which *Leishmania* parasites can escape the host immune response: a signaling point of view. *Clinical Microbiology Reviews* 2005; **18**: 293–305.
4. **Meddeb-Garnaoui A, Zrelli H, Dellagi K.** Effects of tropism and virulence of *Leishmania* parasites on cytokine production by infected human monocytes. *Clinical and Experimental Immunology* 2009; **155**: 199–206.
5. **Herwaldt BL.** Leishmaniasis. *Lancet* 1999; **354**: 1191–1199.
6. **Singh S, Sivakumar R.** Recent advances in the diagnosis of leishmaniasis. *Journal of Postgraduate Medicine* 2003; **49**: 55–60.
7. **Cruz I, et al.** An approach for interlaboratory comparison of conventional and real-time PCR assays for diagnosis of human leishmaniasis. *Experimental Parasitology* 2013; **134**: 281–289.
8. **Reithinger R, Dujardin JC.** Molecular diagnosis of leishmaniasis: current status and future applications. *Journal of Clinical Microbiology* 2007; **45**: 21–25.
9. **Mori Y, Notomi TJ.** Loop-mediated isothermal amplification (LAMP): a rapid, accurate, and cost-effective diagnostic method for infectious diseases. *Journal of Infection and Chemotherapy* 2009; **15**: 62–69.
10. **Karanis P, Ongerth J.** LAMP – a powerful and flexible tool for monitoring microbial pathogens. *Trends in Parasitology* 2009; **25**: 498–499.
11. **Pöschl B, et al.** Comparative diagnosis of malaria infections by microscopy, nested PCR, and LAMP in northern Thailand. *American Journal of Tropical Medicine and Hygiene* 2010; **83**: 729.
12. **Inomata A, et al.** Development and evaluation of a reverse transcription-loop-mediated isothermal amplification assay for rapid and high-sensitive detection of *Cryptosporidium* in water samples. *Water Science and Technology* 2009; **60**: 2167–2172.
13. **Plutzer J, Törökne A, Karanis P.** Combination of ARAD microfibre filtration and LAMP methodology for simple, rapid and cost-effective detection of human pathogenic *Giardia duodenalis* and *Cryptosporidium* spp. in drinking water. *Letters in Applied Microbiology* 2010; **50**: 82–88.
14. **Sotiriadou I, Karanis P.** Evaluation of loop-mediated isothermal amplification for detection of *Toxoplasma gondii* in water samples and comparative findings by polymerase chain reaction and immunofluorescence test (IFT). *Diagnostic Microbiology and Infectious Disease* 2008; **62**: 357–365.
15. **Thekisoe OM, et al.** Detection of *Trypanosoma cruzi* and *T. rangeli* infections from *Rhodnius pallescens* bugs by loop-mediated isothermal amplification (LAMP). *American Journal of Tropical Medicine and Hygiene* 2010; **82**: 855–860.
16. **Xu J, Rong R, et al.** Sensitive and rapid detection of *Schistosoma japonicum* DNA by loop-mediated isothermal amplification (LAMP). *International Journal for Parasitology* 2010; **40**: 327–331.
17. **Tseveleki V, et al.** Cellular FLIP long isoform transgenic mice overcome inherent Th2-biased immune responses to efficiently resolve *Leishmania major* infection. *International Immunology* 2007; **19**: 1183–1189.
18. **Dweik A, et al.** Evaluation of PCR-RFLP (based on ITS-1 and HaeIII) for the detection of *Leishmania* species, using Greek canine isolates and Jordanian clinical material. *Annals of Tropical Medicine and Parasitology* 2007; **101**: 399–407.
19. **Mary C, et al.** Reference values for *Leishmania infantum* parasitemia in different clinical presentations: quantitative polymerase chain reaction for therapeutic monitoring and patient follow-up. *American Journal of Tropical Medicine and Hygiene* 2006; **75**: 858–863.
20. **Notomi T, et al.** Loop-mediated isothermal amplification of DNA. *Nucleic Acids Research* 2000; **28**: E63.
21. **Grimaldi Jr. G, Tesh RB.** Leishmaniasis of the New World: current concepts and implications for future research. *Clinical Microbiology Reviews* 1993; **6**: 230–250.
22. **Takagi H, et al.** Sensitive, specific, and rapid detection of *Leishmania donovani* DNA by loop-mediated isothermal amplification. *American Journal of Tropical Medicine and Hygiene* 2009; **81**: 578–582.
23. **Adams ER, et al.** Development of a reverse transcriptase loop-mediated isothermal amplification (LAMP) assay for the sensitive detection of *Leishmania* parasites in clinical samples. *American Journal of Tropical Medicine and Hygiene* 2010; **82**: 591–596.
24. **Kato H, et al.** Molecular epidemiology for vector research on leishmaniasis. *International Journal of Environmental Research and Public Health* 2010; **7**: 814–826.